

NEW

Introducing the Rifton Adaptive Desk

This desk is the inclusive learning solution that adjusts to every student's needs in seconds. Empower every student with the truly adaptive desk that makes learning accessible for all.



Adaptive Desk

Introducing Rifton's Adaptive Desk – This desk is the inclusive learning solution that adjusts to every student's needs in seconds. With a wide, wheelchair-friendly frame and tool-free continuous height adjustment, it accommodates all learners with ease. The optional 0°–50° surface tilt provides customization for learning activities in seated and standing positions, while the height adjustment makes transitions a breeze. Empower every student with the truly adaptive desk that makes learning accessible for all.



Find important details on the Adaptive Desk at:
www.rifton.com/desk

"We are thrilled to finally have this versatile desk that truly is 'one size fits all.' It can accommodate a wide range of needs, making it an invaluable asset for inclusive learning environments. Whether a child is sitting, standing, using a wheelchair, or a standing system, the desk can be effortlessly adjusted in seconds to provide the optimal support for their unique postural requirements. It's a game-changer that enhances both physical well-being and educational outcomes, making learning more accessible and enjoyable for every child. We couldn't be more impressed with its thoughtful design and the positive impact it has had in our classroom. *"*

Bill Lynch, DPT, Orange-Ulster BOCES, Goshen, NY



Adaptive Desk features

The inclusive learning solution
that quickly adjusts to every student's needs.



Height adjustment:
23" – 43"

Durable, cleanable
desk surface



Optional tilt:
0° – 50°

Foot-activated
continuous height
adjustment

Wide base
accommodates large
Rifton Activity Chair or
powered wheelchair

Therapeutic benefits

1 Sit-to-stand height adjustment

offers variability in learning position for improved posture and physical activity outcomes. Angle adjustment enables ergonomic and visual access to materials and devices.

2 Adjustable elbow supports

provide customized forearm support more effectively than desk cutouts.

3 The mobile base adapts to dynamic classroom environments.



Versatile height and tilt adjustments enable comfortable access in any position.



The desk accommodates most wheelchairs and the large Rifton Activity Chair.





Optional tilt adjustment 0° – 50°



Foot activated continuous height adjustment 23" – 43"



The mobile base adapts to dynamic classroom environments.



The side storage box offers convenient accessory organization.



Adjustable elbow supports offer better forearm support than desk cutouts.




Anti-tip rollers keep the desk from tipping forward.


The Adaptive Desk comes in your choice of surface colors.



Maple

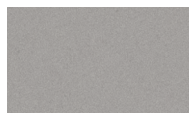

Gray



Black

Components



Tilt

The tilt option allows for continuous adjustment from 0°–50° for ergonomic positioning.



Raised edge

The raised edge keeps book and tablets from sliding off the desk.



Elbow supports

The adjustable elbow supports offer customizable forearm support.



Slide-on ledge

The slide-on ledge offers additional security to keep objects on the desk when tilted.



Device mount clamp adapter

The device mount clamp adapter allows for securing communication devices on raised edge desks.





Casters

The casters provide mobility and can be locked for stability.



Stationary feet

The stationary feet provide a stable base and are the same height as the casters.



Side storage box

The side storage box allows for convenient accessory storage.



Anti-tip rollers

The anti-tip rollers prevent the desk from flipping forward.



Adaptive Desk dimensions

Product dimensions (inches)	R250 Adaptive Desk
Height	23–43
Frame length	36½
Frame depth	28
Top surface dimensions	24 x 40
Weight	68 lb
Max Load	110 lb (55 lb when tilted)
Elbow supports	4½ x 8½, 25 lb limit
Distance between elbow supports	11–20

2025 order form

Scan for sizing and pricing:
<https://www.rifton.com/desk>



		R250
*Frame		R251 \$1,200 <input type="checkbox"/>
*Base	CHOOSE	Casters R252 \$200 <input type="checkbox"/>
		Stationary Feet R253 \$60 <input type="checkbox"/>
Adjustable tilt (requires raised edge)		R254 \$300 <input type="checkbox"/>
*Surface	CHOOSE	Maple with raised edge R261 \$550 <input type="checkbox"/>
		Gray with raised edge R262 \$550 <input type="checkbox"/>
		Black with raised edge R263 \$550 <input type="checkbox"/>
		Maple with plain edge R264 \$400 <input type="checkbox"/>
		Gray with plain edge R265 \$400 <input type="checkbox"/>
Anti-tip rollers		R255 \$120 <input type="checkbox"/>
Side storage box		R256 \$120 <input type="checkbox"/>
Elbow supports		R257 \$120 <input type="checkbox"/>
Device mount adapter (for raised edge)		R258 \$20 <input type="checkbox"/>
Slide-on ledge (requires raised edge)		R259 \$40 <input type="checkbox"/>

* Required